

dermagel[®] micro

The instruction below should be used in conjunction with detailed information on the packaging.

Short description of the product

Sterile, natural rubber latex surgical and protective gloves powder-free, for disposable use

Full description of the product

Reference number	: RC10008060-90
Raw material	: natural rubber latex
Donning powder	: -
External surface	: textured, chlorinated
Internal surface	: polymerized
Cuff	: beaded
Colour	: light brown
Shape	: anatomic, hand specific
Size range	: 6.0; 6.5; 7.0; 7.5; 8.0; 8.5; 9.0
AQL	: 0.65
Sterilization	: Gamma (R)
Packaging	: 1 pair per pouch, 50 pairs per dispenser, 400 pairs per carton
Shelf life	: 3 years (from the date of manufacturing)

Storage instructions

It is recommended to store the gloves in dry place, in the temperature of 5-35°C and to protect them against direct sunlight and fluorescent light.

It is recommended to protect the gloves against humidity. Keep the gloves in a distance of not less than 1m from heating devices, sources of fire and ozone.

MDD classification & compliance

Gloves classified as medical device – class IIa under Council Directive 93/42/EEC.

Compliance with standards: EN ISO 13485:2016, EN ISO 14791:2012, EN ISO 15223-1:2016, EN 1041:2008, EN 455-1:2000, EN 455-2:2009+A2:2013, EN 455-3:2006, EN 455-4:2009, EN 556-1:2001|AC:2006, EN ISO 11737-1:2006, EN ISO 11737-2:2009, EN ISO 11137-1:2015, EN ISO 11137-2:2015, EN ISO 10993-1:2009|AC:2010, EN ISO 10993-5:2009, EN ISO 10993-10:2010.

DNV GL Presafe AS, Notified Body Number 2460:

CE 2460

MDD 93/42/EEC
DNV GL Presafe AS
Veritasveien 3
N-1363 Høyvik, Norway

PPER classification & compliance

Gloves are category III Personal Protective Equipment as per Annex I of the Regulation 2016/425.

Compliance with standards: EN 420:2003+A1:2009, EN ISO 374-1:2016+A1:2018 (type C), EN ISO 374-2:2019, EN 16523-1:2015, EN ISO 374-4:2019, EN ISO 374-5:2016.

EU Type Examination (Module B) and on-going conformity (Module D) Notified Body:

0598

PPER 2016/425
SGS FIMKO OY
Takomotie 8
00380 Helsinki, Finland

Declaration of Conformity can be found under below web address:
<https://en.mercatormedical.eu/>

Intended use

Sterile, surgical and protective gloves, made from natural rubber latex, with thumb positioned towards the palm side of the index finger which reduces the fatigue on the hands, intended to be worn on the hands usually in surgical settings to provide barrier against potentially infectious fluids or other contaminants and to protect the patient by ensuring sterility of the wound environment. These gloves are intended for single use only. Gloves are classified as Medical Devices Class I and as a Personal Protective Equipment Category III. Their design and labelling corresponds to the requirements of the European Regulation 2017/745 on Medical Device and the European Regulation 2016/425 on Personal Protective Equipment. Gloves should be used solely according to their intended application.

Precautions and indications for use

Attention should be drawn to the fact that the confirmed protection level does not reflect the real working conditions where other factors may influence the protection level, for example, abrasion, degradation, atmospheric conditions. Users should take care when using the gloves. Using them solely according to their intended application. These are sterile gloves for single use only and are not intended for cleaning or re-sterilization. The protection is limited to the hand only, the results relate to the palm of gloves and there were tested under laboratory condition.

The gloves should be put on dry hands. Before usage, inspect the gloves for any defect or imperfections. Do not let chemical substances get under the gloves through the cuff. If a chemical substance reaches the skin, wash it away immediately with plenty of water with soap. If the gloves get punctured, torn or broken during their use, take them off and put on the new ones. Avoid using gloves dirty in the inside as they may cause irritation leading to skin inflammation or more serious damages. The gloves should not be used in contact with open fire and to protect against any sharp tools.

The gloves are not intended for welding, electric shock protection, ionizing radiation or from the effect of hot or cold objects.

This information does not reflect the actual duration of protection in the workplace and the differentiation between mixtures and pure chemicals.

The chemical resistance has been assessed under laboratory conditions from samples taken from the palm only (except in case where glove is equal to or over 400 mm – where the cuff is tested also) and relates only to the chemical tested. It can be different if the chemical is used in a mixture.

It is recommended to check that the gloves are suitable for the intended use because the conditions at the workplace may differ from the type test depending on the temperature, abrasion and degradation.

When used, protective gloves may provide less resistance to the dangerous chemical due to changes in physical properties. Movements, snagging, rubbing, degradation caused by the chemical contact etc. may reduce the actual use time significantly. For corrosive chemicals, degradation can be the most important factor to consider in selection of chemical resistant gloves.

Components / hazardous components

Some gloves may contain components known to be a possible cause of allergy for person allergic to them, who may develop contact irritation and/or allergic reaction. Natural rubber latex gloves may cause allergic reactions including anaphylactic reactions. In case of an allergic reaction, seek medical assistance immediately.

Disposal

Used gloves can be contaminated with contagious or other hazardous substances. They should be disposed of in accordance with local regulation. Gloves should be buried or burned under controlled conditions.

Manufacturer

MERCATOR MEDICAL S.A.
ul. H. Modrzejewskiej 30
31-327 Cracow, Poland
www.mercatormedical.eu




Permeation Performance Levels as per EN ISO 374-1:2016+A1:2018					
• Level 1 > 10 min • Level 2 > 30 min • Level 3 > 60 min • Level 4 > 120 min • Level 5 > 240 min • Level 6 > 480 min					
Test results acc. to EN 16523-1:2015		EN ISO 374-4:2019 Degradation [%]	Test results acc. to EN 16523-1:2015		EN ISO 374-4:2019 Degradation [%]
Chemical	Level		Chemical	Level	
40 % Sodium hydroxide	Level 6	-15.6	30% Hydrogen peroxide	Level 3	11.1
99% Isopropanol	Level 0	29.8	37% Formaldehyde	Level 1	20.9
Test acc. to EN ISO 374-2:2019 – Level 3 (ISO 2859)			Test acc. EN ISO 374-5:2016		
Performance Level	AQL		Protection against bacteria & fungi		Pass
Level 3	< 0.65		Protection against viruses		Pass
Level 2	< 1.5				
Level 1	< 4.0				


Symbols which can be used on the packagings

Do not re-use / gloves are intended for single use	Powdered gloves	R Glove is fitting the right hand
Sterilized with Ethylene oxide	Powder free gloves	L Glove is fitting the left hand
Sterilization with Gamma	Presence of polymer coating on the inner surface of the glove	1x para/pair/napa the number of pairs of gloves in the envelope (1)
Do not re-sterilize	Marking of gloves protecting against bacteria and fungi	25 par / pairs / napa the number of pairs of gloves in the dispenser (25)
Raw material – natural rubber latex	Marking of gloves protecting against viruses, bacteria and fungi	50 par / pairs / napa the number of pairs of gloves in the dispenser (50)
Raw material – neoprene	Marking of type A chemical resistant gloves. Six tested chemicals shall be identified by their code letter under pictogram	100 par / pairs / napa the number of pairs of gloves in the carton (100)
Raw material – nitrile	Marking of type B chemical resistant gloves. Three tested chemicals shall be identified by their code letter under pictogram	300 par / pairs / napa the number of pairs of gloves in the carton (300)
Lot / batch number	Marking of type C chemical resistant gloves	400 par / pairs / napa the number of pairs of gloves in the carton (400)
Catalogue number	Keep away from moisture, store in a dry place	1x25 komplet / set / zaimpakovanje the number of sets of gloves in the dispenser (25)
Manufacturer, symbol should be accompanied by name and address of Manufacturer	Keep away from solar and fluorescent light	1x200 komplet / set / zaimpakovanje the number of sets of gloves in the carton (200)
EU Authorised Representative, symbol should be accompanied by name and address of Authorised Representative	Temperature limitation / gloves store in temperature 5-35°C	Consult instructions for use
data produkcji	Keep away from ozone	Do not use, if package is damaged
Expiry date	Package made from paper, qualify for recycling	Package is treated as municipal waste
Indicates compliance with the requirements of Russian market	C PAP	Indicates compliance with the requirements of Ukrainian market

Donning surgical gloves guideline

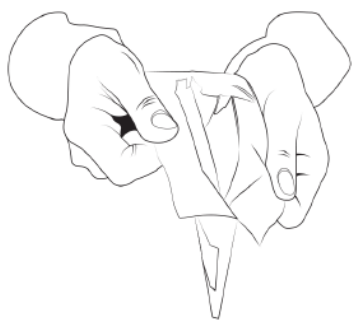

1) Open the packaging and remove the inner sachet containing the gloves.

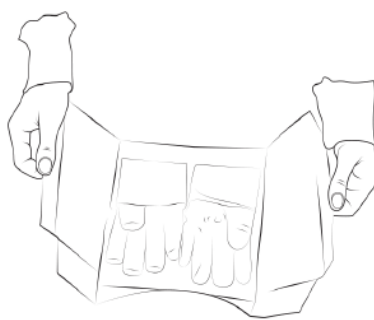

2) Holding it by the wings, open the sachet.



3) Holding the exposed inside cuff of the left glove with the right hand, put on the glove.



4) Insert the left, gloved hand into the cuff opening of the right glove.



5) Pull the glove onto the right hand.



6) Pull the gloves up over the cuffs of the gown.